



## Special Protection

Tracheal compresses can be used, for example, to provide padding between the tube and the stoma; they also absorb secretions coming from the stoma. Other products, such as special cleaning sets and brushes, help to maintain the extremely important high standards of hygiene. A well-cared tracheostoma ensures that the patient feels comfortable and reduces the complications.

Bibs, scarves, and roll-necks can help to replace the missing protective and filtration functions of the nose. They protect the airway from dryness, cold and dust but also act as coverings to improve the patient's appearance.



**REF 940**

**TRACOE care Cleaning Brushes flexible**

Packaging unit: ◆ 10 units

Ordering Number	Bristle-diameter mm	for tube Inside Diameter mm	approx. Total Length mm
940-A	7	5/6	315
940-B	9	7/8	335
940-C	11	9/10	360
940-D	13	11/12	360

Extra-long, specially made for cleaning flexible tracheostomy tubes. The wire used has a rounded tip and the bristles are made of soft nylon.



**REF 915**

**TRACOE care Shower Guard**

Silicone shower protector, contains no DEHP or latex. Soft, adjustable, and can be adapted to different-sized necks.

**Caution: May not be used for full-body bathing or swimming!**



**TRACOE care Tracheostoma Protection**

**REF 908**

Dimensions: ◆ approx. 6 x 6.5 cm

Packaging unit: ◆ 10 units

**REF 909**

Dimensions: ◆ approx. 7.2 x 7.2 cm

Packaging unit: ◆ 10 units

Prevents foreign objects from entering the tracheostoma. With adhesive strip, for single use, for adults and children.



**REF 905**

**TRACOE care Tracheal Compresses**  
aluminium coated, with zigzag slit

Dimensions:           ◆ approx. 9 x 9.8 cm  
Packaging unit:       ◆ 10 units

Tracheal compresses, aluminium coated, avoid the compress sticking to the wound or the skin. They absorb secretions between the stoma and the tube.



**TRACOE care Tracheal Compresses two-ply**

**REF 907-A**   without slit

**REF 907-B**   with zigzag slit

Dimensions:           ◆ approx. 9 x 9.8 cm  
Packaging unit:       ◆ 10 units

Two-ply tracheal compresses absorb a large amount of secretions from the tracheostoma quickly to avoid skin irritations. At the same time they are used as a pad between the tracheostoma and the tracheostomy tube.



**TRACOE care Tracheal Compresses**

**REF 910-A**   without slit

**REF 910-B**   with zigzag slit

Dimensions:           ◆ approx. 9 x 9.8 cm  
Packaging unit:       ◆ 10 units

These tracheal compresses absorb secretions from the tracheostoma to avoid skin irritations.

**REF 957**

**TRACOE care Coated Tracheal Compresses**  
two-ply, with zigzag slit

Dimensions:           ◆ approx. 9 x 10 cm  
Packaging unit:       ◆ 10 units

Two-ply tracheal compresses absorb a large amount of secretions from the tracheostoma quickly to avoid skin irritations. The absorbent coating further avoids the compress sticking to the wound or skin. At the same time they are used as a pad between the tracheostoma and the tracheostomy tube.



**REF 960**

**TRACOE care Coated Tracheal Compresses**  
with zigzag slit

Dimensions:           ◆ approx. 8 x 10 cm  
Packaging unit:       ◆ 10 units

To absorb secretions from the tracheostoma to avoid skin irritations. The coating further avoids the absorbant compress sticking to the wound or the skin.



**REF 961**



**TRACOE care Tracheal Compresses**  
with zigzag slit, for neonates, infants and children

Dimensions:           ◆ approx. 6.5 x 6.3cm  
Packaging unit:       ◆ 10 units

To absorb secretions from the tracheostoma. At the same time they are used as a pad and a protection for the sensitive skin around the tracheostoma. The compresses are ideally adapted to the anatomy of neonates and children.



**REF 969**

**TRACOE care softpad Tracheostomy Stoma Pads oval**



- Dimensions: ♦ approx. 5.8 cm x 3.6 cm  
 Packaging unit: ♦ 20 units  
 ♦ sterile, individually packaged

Oval tracheostomy stoma pads made of open-pored, absorbent polyurethane foam. Avoids decubitus, compensates unevenness due to scars, seals air leaks.

**TRACOE purofoam Tracheostomy Dressings with zigzag slit**



**REF 958** large

- Dimensions: ♦ approx. 9 x 10 x 0.5 cm  
 Packaging unit: ♦ 10 units  
 ♦ sterile, individually packaged

For tracheostomy tube sizes 06 – 12.




**REF 959** small

- Dimensions: ♦ approx. 6.5 x 6.5 x 0.5 cm  
 Packaging unit: ♦ 10 units  
 ♦ sterile, individually packaged

For tracheostomy tube sizes 03 – 05.



**Features of TRACOE purofoam Tracheostomy Dressing REF 958 and REF 959:**

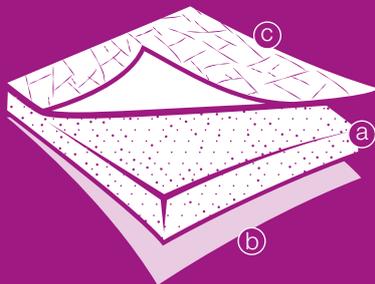

The TRACOE purofoam tracheostomy dressings are made of an absorbent open pore polyurethane foam (a).

The surface layer of the compress that faces the skin allows the absorption of liquids into the foam's core (b). A net-like film stabilises the purofoam dressings and thus reduces the deformation with strong secretion after the moisture absorption.

The opposite side is covered with a beige coating which will prevent leakage of secretions and bacteria through the front of the dressing (c).

Due to the circular opening and the zigzag slit the dressing can be easily replaced with a new one whilst wearing the tracheostomy tube. The TRACOE purofoam tracheostomy dressing can be used to dress a newly established or an existing, irritated tracheostoma, as well as healed but weeping tracheostomas and any other kind of artificial orifice. This dressing meets modern standards of moist wound treatment.